




Quality - Value - Design


# LIFTER SLINGS

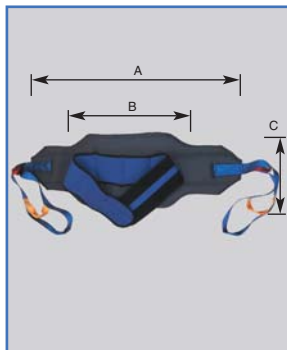

## TOILETING/STANDING SLING

- Suitable for patients who generally have good trunk and hip control
- Easier access for clothing removal during toileting and dressing
- Allows for safe comfortable chair and bed transfers
- Elasticated neoprene waist strap
- Allows patients to be easily brought to a standing position
- Available in Polyester-Nylon only
- Positioning straps
- Colour coded attachments straps
- N.A.T.A.\* Tested safe working load of 250Kg
- Made to Australian Standards
- Australian made and owned

	A	B	C	Poly
Small/Medium	1250	400	800	STSP1
Large/X-large	1300	500	900	STSP2
Special				STSP9

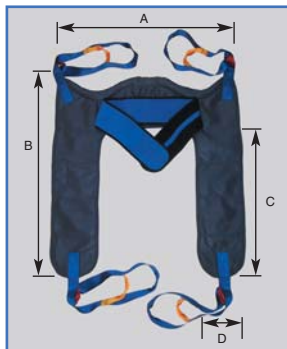

## COMMUNE SLING

- Elasticated neoprene waist strap
- Long divided leg pieces
- Easy access for clothing removal during toileting and dressing
- Soft padding throughout sling
- Suitable for patients who have generally good trunk and hip control
- Available in Polyester – Nylon only
- Positioning straps
- Colour coded attachment straps
- N.A.T.A.\* Tested safe working load of 250Kg
- Made to Australian standards
- Australian made and owned


	A	B	C	D	Poly
Small	800	900	650	200	SCSP1
Medium	850	900	650	200	SCSP2
Large	950	900	650	200	SCSP3
Special					SCSP9



## GENERAL PURPOSE SLING

- Suitable for patients with good head control
- Quick and easy to fit
- Divided legs can be crossed in the centre, uncrossed or positioned under thighs
- Available in Polyester-Nylon and Polyester-Mesh
- Positioning straps
- Colour coded attachment straps
- N.A.T.A.\* Tested safe working load of 250Kg
- Made to Australian standards
- Australian made and owned

	A	B	C	D	Poly	Mesh
Child	750	900	500	300	SGPP0	SGPM0
Small	800	1150	550	400	SGPP1	SGPM1
Medium	900	1300	600	450	SGPP2	SGPM2
Large	950	1350	650	500	SGPP3	SGPM3
Special					SGPP9	SGPM9

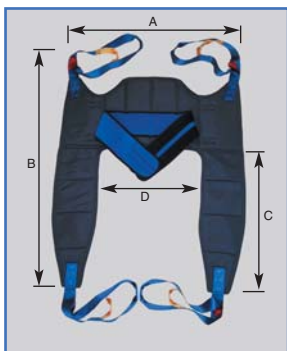

## DELUXE COMMUNE SLING

- Soft padding throughout the sling
- Excellent back support and comfort
- Suitable for patients who generally have good trunk and hip control
- Long divided legs
- Elasticated neoprene waist strap
- Easy access for clothing removal during toileting
- Available in Polyester-Nylon only
- Positioning straps
- Colour coded attachment straps
- N.A.T.A.\* Tested safe working load of 250Kg
- Made to Australian standards
- Australian made and owned


	A	B	C	D	Poly
Small	750	950	600	400	SCDP1
Medium	850	1000	650	500	SCDP2
Large	950	1010	750	600	SCDP3
Special					SCDP9

# care and washing

CARE AND WASHING INSTRUCTIONS OF SLINGS

- Do not lift a patient under any circumstances if there are any sling defects
- The Velcro straps must remain closed when not in use
- Periodically inspect the slings for defects
- Repair any sign of tears immediately
- Do not use positioning straps to lift patient, it is only to guide patient
- Washing the slings at a temperature greater than 75° Degrees Celcius will reduce the life of the sling and may cause shrinkage or deformation of the foam
- Dry thoroughly before using sling and dry in away from direct heat and sunlight
- Washing instructions are placed on each sling and should be adhered to at all times

## STANDARD STRAP DIMENSIONS

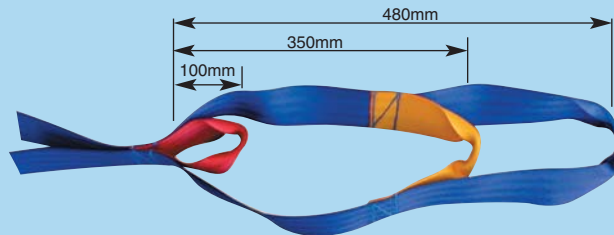

## TITAN SOLID

Description:	Sling / Medical Fabric
Yarn:	100% Polyester
Yarn Type:	167 D/Tex Flat Polyester
Yarn Tenacity:	Medium
Weight:	Minimum 300 gms
Width:	127 / 132 cm
Bursting Strength:	Minimum 78 p.s.i. (7" Mullins)
Certificate:	Yes to Uk/SC/5188A Patter

## TITAN MESH

Description:	Sling / Medical Fabric
Yarn:	100% Polyester
Yarn Tenacity:	High Tenacity Bright
Weight:	Minimum 375 gms
Width:	127 / 132 cm
Bursting Strength:	Warp 1000N                      BS 2576 Weft 1000N                      BS 2576
Other:	Warp Knitted Wales P / Inch Minimum 17 Courses P / Inch              Minimum 23
Certificate:	Yes to UK/SC/5188A Pattern 8795B
Recommended Washing:	Standard disinfection Wash Cycle up to 75 Degrees Celsius

All the above values are of an average and are subject to a certain variation in bulk production  
Titan Solid & Titan Mesh material is available in Black & Navy

\* National Association of Testing Authorities Australia (NATA)

WARRANTY: 2 years. Conditions apply

We reserve the right to alter specifications or finishes without prior notice in accordance with our policy to improve product design.

For more information contact:

**Jomor Healthcare P/L** 308 Lower Dandenong Road  
Mordialloc, VIC 3195  
**ph:** 1300 651 235 **fax:** (03) 9588 2255  
**email:** sales@jomor.com.au **web:** www.jomor.com.au

**JOMOR**  
**HEALTHCARE**  
*Quality - Value - Design*

